




It goes everywhere, does it all, and thinks the way you do!



- **Best argon beam.** Robust arc length without sacrificing limited tissue depth.
- **Truly standardizable.** The only generator and pump you need. Does it all: argon coagulation, pulse cut for sphincterotomy, true GI bipolar, TouchSoft® coagulation and tried and tested polypectomy, EMR, and ESD outputs.
- **Fits everywhere.** Patented, compact argon canister system goes inside the all in one unit. Fits on any boom, travel or endo cart. Put it everywhere you need it!
- **Disposable argon gas canister** eliminates bulky argon! Compact, easy, economical.
- **Lavage included.** The only lavage pump you need!
- **A revolution in ease of use!** It thinks the way you do! Touch, confirm and treat! Simpler to be safer. Foot pedal is either power or water. No more color guessing!
- **Tap a Tool®** Just a touch gives you default confidence for endo specific tools.
- **Designed just for endoscopy.** No confusing adapters or ports or programs!
- **Most advanced tissue sensing:** Auto output control for constancy in tissue results.
- **PSS--** Genii exclusive **Pad Safe System** for advanced grounding pad safety with revolutionary three color pad sense and no light cue for bipolar.
- **Standby** for greater safety.
- **Affordable:** you can afford one for every room, every cart. You CAN now achieve the cost benefits of being truly standardized! **It's not just smaller--it's smarter!**

The revolutionary *gi4000* has four output methods and seven monopolar modes, for perfect confidence in all procedures—basic to the most complex.

Methods: Argon Coagulation, Monopolar, Bipolar, and Lavage

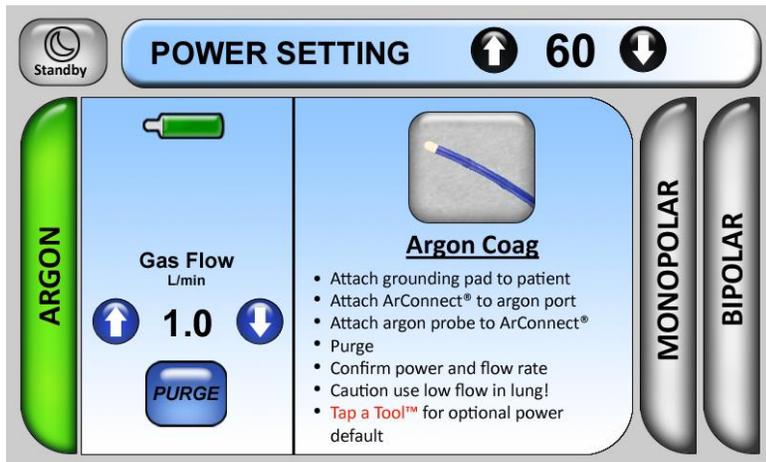

Argon:

- Proprietary ArC Smart™ beam for robust arc length without sacrificing limited tissue depth, and special voltage control to reduce video interference
- Revolutionary disposable argon canister
- Quick, reliable arc
- Easy set up with disposable ArConnect®. Probe choice without reprocessing




Seven Monopolar Outputs

- **TouchSoft®** low voltage coagulation for gentle, depth limited coagulation.
- **Coag:** tried and trusted for standard snare polypectomy.
- **Blend Coag** for more cut, larger polyps, EMR
- **Pulse Blend Cut:** controlled cutting with adequate hemostasis when desired for sphincterotomy, EMR, ESD
- **Blend Cut:** like pulse blend but not fractionated (no pulse)
- **Pulse Cut:** controlled pulsed cutting with minimal coag for sphincterotomy, EMR, ESD
- **Cut:** constant maximum cut, minimum coagulation

Bipolar

- Best endo bipolar

Lavage—always ready!

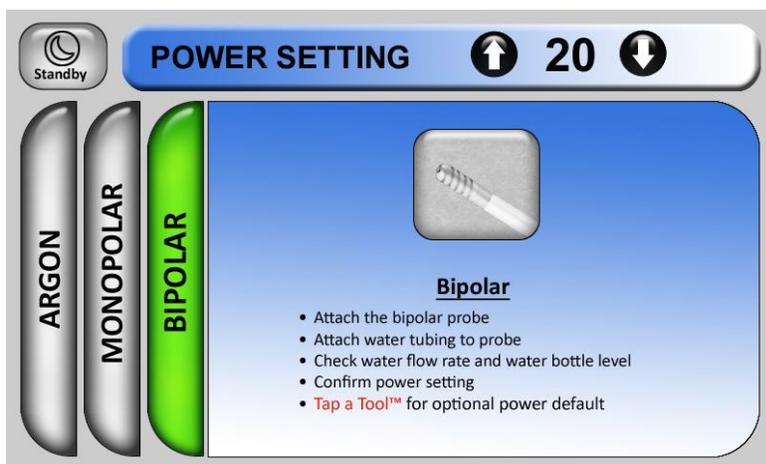

Find out more at www.genii-gi.com or 855-501-4810

The *gi4000* is made in the USA

Rev 5 1 15